



Method number: 107

Matrix: Air

	Maneb	Zineb
Target concentration:	5 mg/m ³ TWA	5 mg/m ³ TWA
OSHA PEL:	None	None
ACGIH TLV:	None*	None

* 5 mg/m³ for manganese dust and compounds, as Mn (1994-95)
0.2 mg/m³ for elemental manganese and inorganic compounds (1995-96)

Procedure: Samples are collected by drawing known volumes of air through sampling cassettes containing a membrane filter made of mixed esters of cellulose. Samples are extracted with an aqueous solution of 5% cysteine and 5% EDTA and analyzed by LC using a UV detector.

Recommended air volume and sampling rate: 500 L at 2.0 L/min

	Maneb	Zineb
Reliable quantitation limit:	220 µg/m ³	103 µg/m ³
Standard error of estimate:	10.1%	9.9 %

Status of method: Evaluated method. This method has been subjected to the established evaluation procedures of the Organic Methods Evaluation Branch.

Date: February 1996

Chemist: Yihlin Chan


Organic Methods Evaluation Branch
OSHA Salt Lake Technical Center
Salt Lake City, UT 84165-0200

1. General Discussion

1.1 Background

1.1.1 History

Thiocarbamates and dithiocarbamates are two of the most important families of fungicides (Ref. 5.1). The latter family includes compounds such as ziram, nabam, ferbam, maneb, and zineb. Of these, maneb and zineb are the most commonly used. They and nabam are ethylenebisdithiocarbamates (EBDTC). Although their toxicities are not very high (oral LD₅₀s for mouse are 4100 mg/kg and 7600 mg/kg for maneb and zineb, respectively), monitoring of EBDTC in the working environment or in food residues is important because they can decompose on heating to form ethylene thiourea, an animal carcinogen (Ref. 5.2).

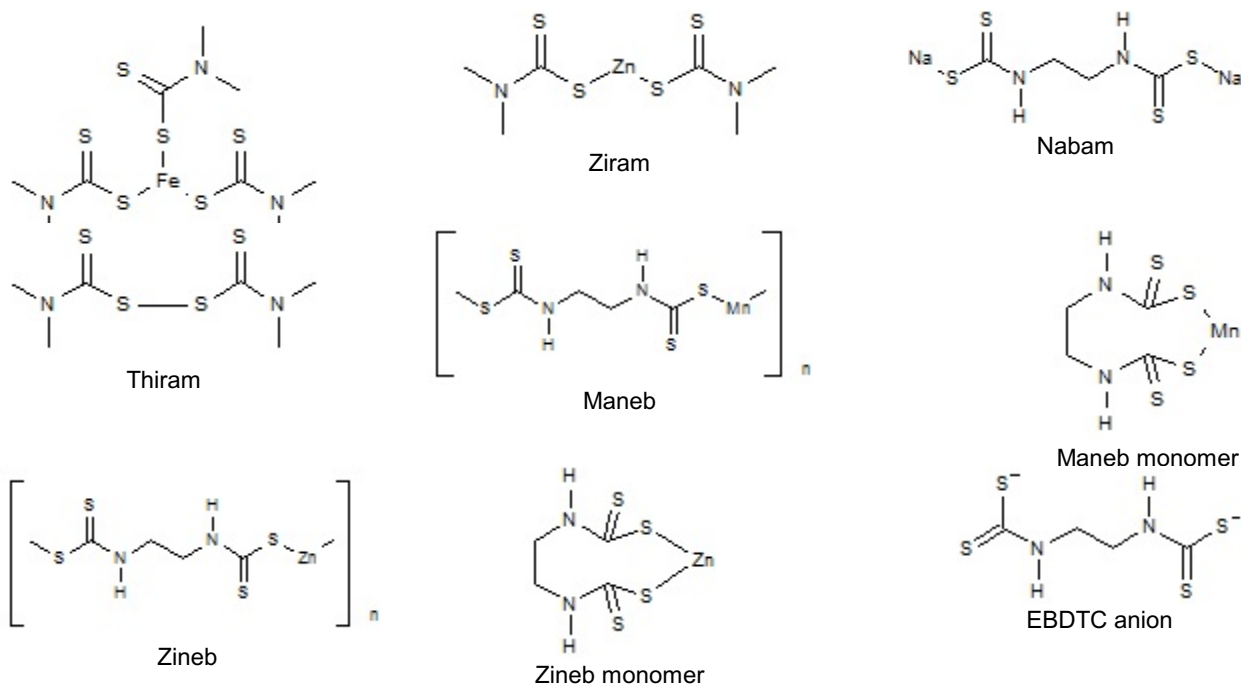

Currently the most common analytical method for maneb and zineb is based on acid hydrolysis with the measurement of the released carbon disulfide by head-space GC analysis (Ref. 5.3). The method is nonspecific because it does not distinguish dialkyl dithiocarbamates (such as ziram and ferbam) or thiuram disulfide (such as thiram) from EBDTC.

When OSHA SLTC first received a set of samples for maneb analysis, it was noted that maneb is soluble in chloroform according to Merck Index. The samples were extracted with chloroform and analyzed by normal-phase LC. The maneb in the eluted peak was confirmed by mass spectrometry (direct probe). Later it was decided to do a fully-validated method for maneb and zineb, but a more detailed literature search revealed it was not so simple.

Although Merck Index listed their structures as monomeric, maneb and zineb are really polymeric. This makes them difficult to analyze because they are not soluble in most solvents. The reason Merck Index (Ref. 5.4) listed their structures as monomeric and described them as soluble in chloroform and other solvents was probably because the scientists who first synthesized these compounds purified them by recrystallization and consequently obtained monomeric maneb and zineb. In the eleventh edition of the Merck Index, the structures of maneb and zineb were corrected to the polymeric forms but they are still listed as soluble in chloroform etc. (Ref. 5.5)

One way to dissolve maneb or zineb is to use chelating agents such as ethylenediamine-tetraacetic acid (EDTA) to strip away the bivalent metals (zinc or manganese) and release the EBDTC anion. The released EBDTC is in the form of sodium salt, namely, nabam. Nabam has been directly analyzed by ion chromatography using UV detector at 286 nm (Ref. 5.6). Miles and Zhou acid-hydrolyzed nabam to ethylenediamine and fluorogenically labeled the latter with *o*-phthalaldehyde-mercaptoethanol (Ref. 5.7). Unfortunately, nabam decomposes rapidly in solution. Bardarov and Zaikov used ascorbic acid as a stabilizer (Ref. 5.6). Kunugi used cysteine as a stabilizer in his analysis of maneb and zineb residues in animal feed (Ref. 5.8). Others have converted the EBDTC anion into methyl ester to avoid the instability problem (Ref. 5.9).

When developing this method, DMSO and DMF were first tested to see if they dissolve and depolymerize maneb and zineb. Maneb and zineb went into solution readily but LC analysis under various conditions were unsuccessful. Nabam was next synthesized from ethylenediamine and carbon disulfide. With the authentic nabam in hand, the LC conditions for its analysis were developed, and nabam was found to decompose when dissolved in DMSO.

In considering the analytical procedure, the head-space GC analysis of carbon disulfide after acid hydrolysis of EBDTC was considered cumbersome. Besides, it suffers from many interferences. Derivatization of the EBDTC anion to methyl ester involves liquid-liquid extraction which is to be avoided if possible. In the end, maneb and zineb were converted to nabam and the EBDTC anion analyzed, even though this meant that the two analytes could not be differentiated. With respect to the stability of EBDTC in solution, the extraction solutions of Bardarov and Kunugi were compared, and Kunugi's formula was found to be much better. The stock solutions prepared in the Bardarov's solvent began to decompose within hours, while those prepared in the Kunugi's stayed unchanged for more than 5 days. This is true even for those stocks prepared in the 49-day-old solvent (Section 4.12).

Because maneb and zineb cannot be dissolved without decomposition, many of the standard OSHA tests for validating a method cannot be followed. The analyte cannot be accurately liquid spiked onto the sampling medium. The aerosol generator could not be used with either an atomizer or a vibrating orifice because these require the source material be dissolved in some kind of solvent. Marple and coworkers developed a dust generator in which the powder is fed by a bead-chain conveyor into a fluidized bed where it is deagglomerated and aerosolized (Ref. 5.10). This kind of dust generator is available from TSI Incorporated of St. Paul, Minnesota (Model 3400). But with our set up for the generator and the test chamber we were unable to obtain a steady, uniform dust atmosphere.


Considering the dusty nature of the analytes, an electrically conductive carbon-filled polypropylene 3-piece sampling cassette with a 51-mm extension, containing a support pad and a 0.8 μm membrane filter made of mixed esters of cellulose was selected for the sampling medium. The test samples were prepared by weighing maneb or zineb directly onto the filter. This would mean that numbers such as the reliable quantitation limit can only be as good as the accuracy of the balance. Mixtures of maneb or zineb with sucrose (10.0% maneb or 11.0% zineb by weight) were used as the test materials in order to extend the reference materials and to improve precision. Sucrose was selected because it was sometimes used as a wetting powder in the field. The mixtures were prepared by using a freezer mill, where maneb (or zineb) and sucrose in a tube are pounded rapidly with a stainless steel ball while submerged in liquid nitrogen.

1.1.2 Toxic effects (This section is for information only and should not be taken as the basis of OSHA policy.)

Maneb can cause irritation of eyes, nose, and throat. Under normal use conditions, maneb is generally regarded as harmless, except for occasional signs of local irritation. However,

it reportedly caused acute renal failure and ECG abnormalities in a 62-year-old man exposed while applying the compound to his garden (Ref. 5.11). Maneb has been tested in mice and rats by oral administration and by single subcutaneous injection. Oral administration produced an increased incidence of lung tumors in mice of one strain, but no increase was observed in three other strains. The studies in rats cannot be evaluated due to the small number of surviving animals. IARC was unable to make an evaluation of the carcinogenicity of mane. (Ref. 5.2)

Zineb can cause irritation to eyes, nose, and throat and is harmful if inhaled. Zineb produced an increased incidence of lung tumors after its oral administration in one strain of mice. Systemic reticulum-cell sarcomas were observed in mice and a variety of sarcomas in rats after its subcutaneous administration. No increases in tumor incidences were observed in two other strains of mice and in two limited studies in rats following oral administration. The available data do not allow an evaluation of the carcinogenicity of zineb to be made. (Ref. 5.2)

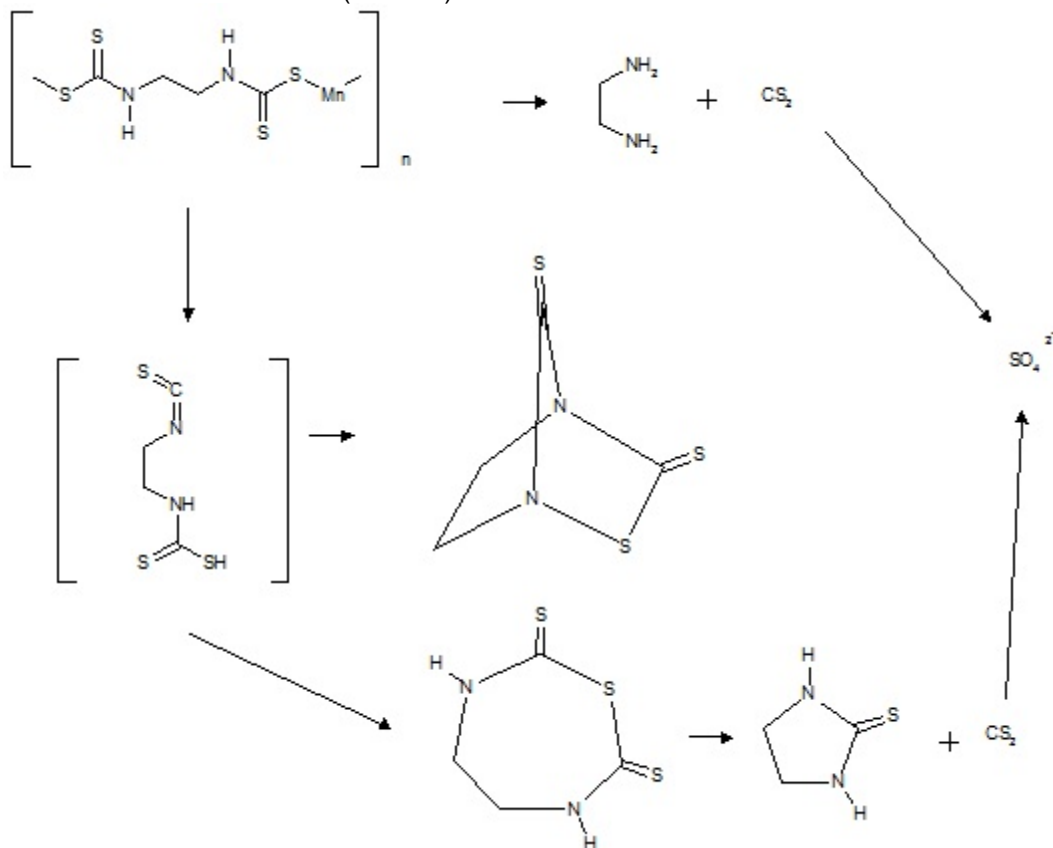

In rats, 55% of an oral dose of mane was excreted in the urine and feces within three days. After 24 hours, the body organs contained 1.2% of the dose as metabolites, and on day 5, less than 0.18%. Ethylenediamine, ethylenebisthiuram monosulfide and ethylene thiourea were present in the urine and feces. Ethylene thiourea is a known animal carcinogen. The above metabolic pathways had been suggested for mane and zineb in rats. (Ref. 5.2).

1.1.3 Workplace exposure

Maneb is used exclusively as a broad spectrum contact fungicide and is registered for use on more than 46 crops in the United States. The principal diseases controlled by mane are early and late blight of potato and tomato, downy mildew and anthracnose on a number of vegetables and the so-called 'rot' diseases of fruits such as apricots, peaches and grapes. It is also used for seed treatment of small grains such as wheat. (Ref. 5.2)

Workers handling various formulations of maneb in the applications mentioned above may be exposed.

Zineb is a fungicide registered for use in the United States on more than 50 crops, including fruits, vegetables, ornamental plants, and for treatment of many seeds. Zineb is also registered for use as a fungicide in paints and for mold control on fabrics, leather, paper, plastic and wood surfaces. (Ref. 5.2) Workers handling zineb in its various formulations in the applications mentioned above may be exposed.


An acceptable daily intake for man of 0-0.005 mg per kilogram body weight for all dithiocarbamate fungicides was established jointly by the Food and Agriculture Organization and the World Health Organization in 1974. (Ref. 5.2)

1.1.4 Physical properties and other descriptive information (Ref. 5.2)

Maneb

CAS no.: 12427-38-2
synonyms: 1,2-ethanediybis(carbamodithioato)(2-)-manganese; 1,2-ethanediybis(carbamodithioic acid, manganese complex); 1,2-ethanediybismaneb, manganese (2+) salt (1:1); 1,2-ethylenediybis(carbamodithioato)manganese; ethylenebis(dithiocarbamic acid), manganese salt; Chem Neb; Chloroble M; CR 3029; Dithane M22; ENT 14875; ethylenebis(dithiocarbamate), manganese; ethylenebis(dithiocarbamic acid) manganous salt; ethylenebis(dithiocarbamate)manganese; F 10; Kypman 80; Lonocol M; Manam; Maneba; Manebgan; Manesan; Manganese (II) ethylenedi(dithiocarbamate); manganese ethylene-1,2-bisdithiocarbamate; Manzate; Nereb; Nesper; Plantifog 160M; Polyram M; Rhodianebe; Sopranebe; Tersan-LSR; Trimangol; Tubothane; Maneb 80; manganous ethylenebis(dithiocarbamate); Trimangol 80; Aamangan; Maneb ZL4; Manzate 200; M-Diphar; MnEBD; MEB; Remasan chloroble M; manganese ethylenebis(dithiocarbamate)

structural formula:

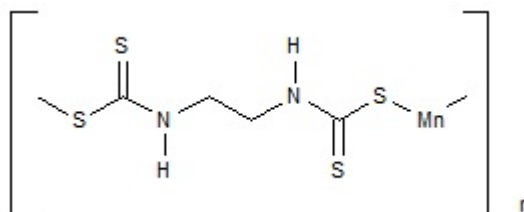


formula wt: 265.29
melting point: decomposes before melting
appearance: yellow-brownish powder
specific gravity: 1.92
vapor pressure: less than 1×10^{-5} Pa at 20°C
solubility: insoluble in most solvents; “soluble in chloroform, pyridine; moderately soluble in water” - (Merck Index. Probably for monomeric maneb.)

Zineb

CAS no.: 12122-67-7
synonyms: 1,2-ethanediybis(carbamodithioato)(2-)-zinc; Aaphytora; Aphytora; Aspor; Asporum; Bercema; Blizene; Carbadine; 1,2-ethanediybis(carbamodithioic acid, zinc salt; Crittox; Daisen;

Deikusol; Discon; Dithane 65; Dithane Z; Dithane Z-78; Ethyl Zimate; ethylenebis(dithiocarbamate)zinc; ethylenebis(dithiocarbamic acid)zinc salt; Fungo-Pulvit; Hexathane; Kypzin; Lipotan; Lirotan; Lonacol; Lonocol; Micide 55; Novozir; Novozir N; Parzate Zineb; Perozin; Perozine; Perozine 75B; Pilzol SZ; Thionic M; Tiezene; Unizeb; Zebenide; Zebtox; zinc ethylenebis(dithiocarbamate); zinc N,N'-ethylenebis(dithiocarbamate); zinc N,N'-ethylenebis(dithiocarbamate); ((1,2-ethylenebis(carbamodithioato))(2-))zinc; Zineb 80; Zinosan

structural formula:

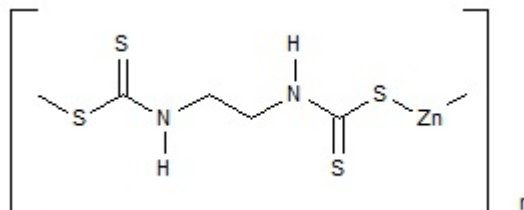

formula wt:	275.75
appearance:	light tan powder
melting point:	decomposes on heating
specific gravity:	1.92
vapor pressure:	negligible at 25°C
solubility:	insoluble in most solvents; "soluble in carbon disulfide, chloroform, pyridine; practically insoluble in water" - (Merck Index, probably for monomeric zineb.)

The analyte air concentrations throughout this method are based on the recommended sampling and analytical parameters.

1.2 Limit defining parameters

1.2.1 Detection limit of the analytical procedure

The detection limits of the analytical procedure are 0.84 and 0.86 ng on column for maneb and zineb, respectively. These are the amounts of analytes that will give responses that are significantly different from the background responses of reagent blanks. (Sections 4.1 and 4.2)

1.2.2 Detection limit of the overall procedure

The detection limits of the overall procedure are 33 and 15 µg per sample (66 and 30 µg/m³) for maneb and zineb, respectively. These are the amounts of analyte spiked on the sampler that will give responses that are significantly different from the background responses of sampler blanks. (Sections 4.1 and 4.3)

1.2.3 Reliable quantitation limit

The reliable quantitation limits are 110 and 51 µg per sample (220 and 103 µg/m³) for maneb and zineb, respectively. These are the amounts of analyte spiked on a sampler that will give signals that are considered the lower limits for precise quantitative measurements. (Section 4.4)

1.2.4 Precision (analytical procedure)

The precision of the analytical procedure, measured as the pooled relative standard deviations over a concentration range equivalent to 0.5 to 2 times the target concentration, are 0.64% and 0.34% for maneb and zineb, respectively. (Section 4.5)

1.2.5 Precision (overall procedure)

The precision of the overall procedure at the 95% confidence level for the ambient temperature 15-day storage tests (at the target concentration) are $\pm 19.8\%$ and $\pm 19.4\%$ for maneb and zineb, respectively (Section 4.6). These include additional 5% for sampling error.

1.2.6 Recovery

The recovery of ethylenebisdithiocarbamate from samples used in 15-day storage tests remained above 93.1% and 96.5% for maneb and zineb, respectively, when the samples were stored at ambient temperature. (Section 4.7)

1.2.7 Reproducibility

Twelve samples, prepared by weighing, were submitted to an SLTC organic service branch for analysis, using a draft copy of this procedure. The samples were analyzed after 2 days of storage at ambient temperature. No individual sample result deviated from its theoretical value by more than the precision reported in Section 1.2.5. (Section 4.8)

2. Sampling Procedure

2.1 Apparatus

2.1.1 A personal sampling pump, calibrated to $\pm 5\%$ of the recommended flow rate with the sampling device attached.

2.1.2 A conductive 3-piece sampling cassette with a 51-mm extension, containing a support pad (25-mm diameter) and a 0.8 μm membrane filter made of mixed esters of cellulose. The sampling media used in this study were obtained from Gelman (catalog number 4375). It contained a GN-4 filter made of mixed esters of cellulose.

2.2 Reagents

None required.

2.3 Technique

2.3.1 Remove the top piece of the cassette for open-face sampling.

2.3.2 Attach the sampler to the sampling pump with a piece of flexible tubing and place it in the worker's breathing zone.

2.3.3 Air should not pass through any hose or tubing before entering the sampling cassette.

2.3.4 After sampling replace the top piece and cap both ends. Wrap each sample with a Form OSHA-21 seal.

2.3.5 Record air volume for each sample.

2.3.6 Submit at least one blank with each set of samples. Blanks should be handled in the same manner as samples, except no air is drawn through them.

2.3.7 List any compounds that could be considered potential interferences.

2.4 Sampler capacity

Sampling capacity was not tested. Generally dusts of such low vapor pressure as maneb or zineb are not expected to show significant loss due to evaporation or sublimation. Retention efficiencies were tested, with the recoveries of 97.0% and 102.8% for maneb and zineb, respectively (Section 4.9).

2.5 Extraction efficiency

2.5.1 The average extraction efficiencies for ethylenebisdithiocarbamate from mixed cellulose ester membrane filters spiked with maneb/sucrose or zineb/sucrose mixture, over the range of 0.5 to 2.0 times the target concentration, were 100.4% and 98.5% for maneb and zineb, respectively. (Section 4.10.1)

2.5.2 The extraction efficiencies at 0.05, 0.1, and 0.2 times the target concentration (TC) are listed below. (Section 4.10.1)

Table 2.5.2
Extraction Efficiencies (%) at 0.05, 0.1, and 0.2 Times the Target Concentration

	maneb	zineb
0.05× TC	96.9	96.5
0.1× TC	100.6	102.7
0.2× TC	92.3	100.6

2.5.3 Extracted samples remain stable for at least 24 h. (Section 4.10.2)

2.6 Recommended air volume and sampling rate

2.6.1 The recommended air volume is 500 L at 2.0 L/min.

2.6.2 For short-term sampling the recommended air volume is 30 L at 2.0 L/min.

2.6.3 When short-term samples are collected, the air concentrations equivalent to the reliable quantitation limits become larger. For example, the reliable quantitation limit is 3.7 mg/m³ for maneb when 30 L is collected.

2.7 Interferences (sampling)

None.

2.8 Safety precautions (sampling)

2.8.1 The sampling equipment should be attached to the worker in such a manner that it will not interfere with work performance or safety.

2.8.2 All safety practices that apply to the work area being sampled should be followed.

3. Analytical Procedure

3.1 Apparatus

3.1.1 An LC equipped with a UV detector. A BAS 200 HPLC (Bioanalytical Systems, Inc., West Lafayette, Indiana) equipped with a UV detector and a Waters 712 autosampler were used in this evaluation.

3.1.2 An anion-exchange column capable of separating ethylenebisdithiocarbamate from any interferences. A Hamilton PRP-X100 column (4.1 mm × 150 mm) column was used in this evaluation.

3.1.3 An electronic integrator or other suitable means of measuring detector response. A Waters 860 Networking Computer System was used in this evaluation.

3.1.4 Glass vials, 20-mL, with poly(tetrafluoroethylene)-lined caps for extracting samples.

3.1.5 A dispenser capable of delivering 4.0 mL of extraction solution.

3.2 Reagents

3.2.1 Maneb. Maneb, 95%, was obtained from Chem Services.

3.2.2 Zineb. Zineb, Tech grade, was obtained from Chem Services.

3.2.3 Ethylenediaminetetraacetic acid, tetrasodium salt hydrate. Ethylenediaminetetraacetic acid, tetrasodium salt hydrate, 98%, was obtained from Aldrich.

3.2.4 *L*-Cysteine hydrochloride hydrate. *L*-Cysteine hydrochloride hydrate, 99%, was obtained from Aldrich.

3.2.5 Sodium perchlorate. Sodium perchlorate, HPLC grade, was obtained from Fisher.

3.2.6 Sodium hydroxide. Sodium hydroxide, reagent grade, was obtained from VWR.

3.2.7 Extraction solution. Dissolve 50 g of ethylenediaminetetraacetic acid tetrasodium salt hydrate and 50 g of *L*-cysteine hydrochloride hydrate in 800 mL of water. Adjust to pH 9.6 with 12 *N* sodium hydroxide. Make the final volume to 1000 mL with water. Store in a brown bottle and use within a month.

3.2.8 LC mobile phase. Dissolve 3.8 g of ethylenediaminetetraacetic acid tetrasodium salt hydrate and 8.4 g of sodium perchlorate in 1000 mL of water.

3.3 Standard preparation

3.3.1 Prepare stock standards by dissolving weighed amounts of mane b or zine b in the extraction solvent and sonicating for 60 min.

3.3.2 Prepare analytical standards by diluting the stock standards with the extraction solvent. For mane b or zine b, a 625 µg/mL standard solution corresponds to the target concentration.

3.3.3 Prepare a sufficient number of analytical standards to generate a calibration curve. Analytical standard concentrations must bracket sample concentrations.

3.4 Sample preparation

3.4.1 Transfer the filter with its collected dust to a glass vial. Discard the supporting pad.

3.4.2 Add 4.0 mL of the extraction solvent to each vial.

3.4.3 Cap the vials and shake them on a mechanical shaker for 60 min.

3.5 Analysis

3.5.1 HPLC conditions

column: Hamilton PRP-X100 (150 mm, 4.1-mm i.d., 10- μ m particle size)
mobile phase: 0.01 M ethylenediaminetetraacetic acid, tetrasodium salt, 0.06 M sodium perchlorate
flow rate: 1.0 mL/min
UV detector: 286 nm
injection size: 10 μ L
retention time: ethylenebisdithiocarbamate 5.0 min

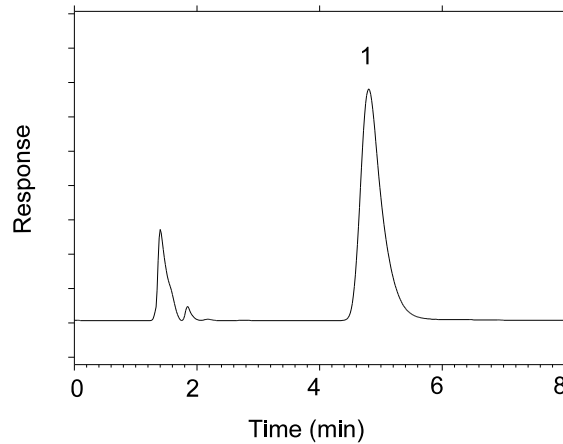


Figure 3.5.1. Chromatogram at target concentration. Key: 1 = ethylenebisdithiocarbamate.

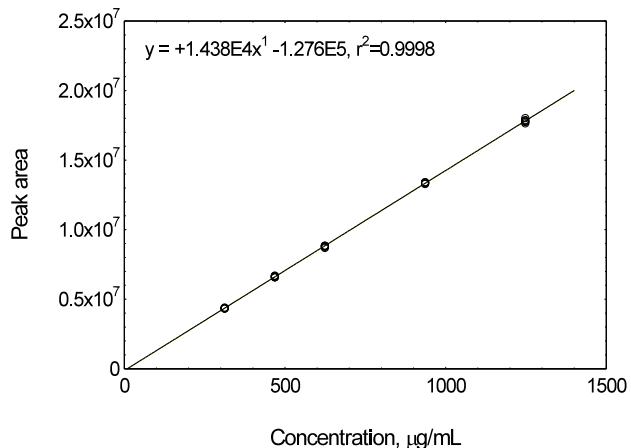


Figure 3.5.3.1. Calibration curve of maneb.

3.5.2 Measure peak areas by an electronic integrator or other suitable means.

3.5.3 Prepare a calibration curve by plotting micrograms per milliliter versus peak areas of standards. Bracket the samples with analytical standards.

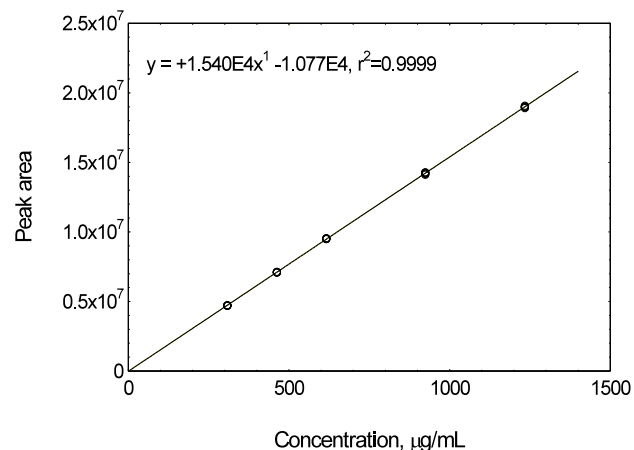


Figure 3.5.3.2. Calibration curve of zineb.

3.6 Interferences (analytical)

3.6.1 Nabam and mancozeb (a coordination complex of zinc ion with manganese ethylenebisdithiocarbamate) cause positive interferences because they produce ethylenebisdithiocarbamate when dissolved in the extraction solvent.

3.6.2 Any other compound that absorbs at 286 nm and has a similar retention time as ethylenebisdithiocarbamate is a potential interference. If any potential interferences were reported, they should be considered before samples are extracted. Generally, chromatographic conditions can be altered to separate an interference from the analyte.

3.6.3 When necessary, the identity or purity of an analyte peak may be confirmed with additional analytical data (Section 4.11).

3.7 Calculations

The amount (in micrograms) of maneb or zineb per milliliter is obtained from the appropriate calibration curve. This amount is corrected by subtracting the amount (if any) found in the blank. The air concentration is calculated using the following formula.

$$\text{mg/m}^3 = \frac{(\text{micrograms per mL}) \times (\text{extraction volume, mL})}{(\text{liters of air sampled}) \times (\text{extraction efficiency})}$$

where: Extraction volume = 4 mL

Extraction efficiency = 1.004 for maneb, or 0.985 for zineb

3.8 Safety precautions (analytical)

3.8.1 Adhere to the rules set down in your Chemical Hygiene Plan.

3.8.2 Avoid skin contact and inhalation of all chemicals.

3.8.3 Wear safety glasses and a lab coat at all times while in the lab area.

3.8.4 Sodium perchlorate is a strong oxidizer.

4. Backup Data

4.1 Determination of detection limits

Detection limits (DL), in general, are defined as the amount (or concentration) of analyte that gives a response (Y_{DL}) that is significantly different (three standard deviations (SD_{BR})) from the background response (Y_{BR}).

$$Y_{DL} - Y_{BR} = 3(SD_{BR})$$

The direct measurement of Y_{BR} and SD_{BR} in chromatographic methods is typically inconvenient and difficult because Y_{BR} is usually extremely low. Estimates of these parameters can be made with data obtained from the analysis of a series of analytical standards or samples whose responses are in the vicinity of the background response. The regression curve obtained for a plot of instrument response versus concentration of analyte will usually be linear. Assuming SD_{BR} and the precision of data about the curve are similar, the standard error of estimate (SEE) for the regression curve can be substituted for SD_{BR} in the above equation. The following calculations derive a formula for DL:

$$SEE = \sqrt{\frac{\sum (Y_{obs} - Y_{est})^2}{n - k}}$$

Y_{obs} = observed response
 Y_{est} = estimated response from regression curve
 n = total no. of data points
 $k = 2$ for a linear regression curve

At point Y_{DL} on the regression curve

$$Y_{DL} = A(DL) + Y_{BR} \quad A = \text{analytical sensitivity (slope)}$$

therefore

$$DL = \frac{(Y_{DL} - Y_{BR})}{A}$$

Substituting $3(SEE) + Y_{BR}$ for Y_{DL} gives

$$DL = \frac{3(SEE)}{A}$$

4.2 Detection limit of the analytical procedure (DLAP)

The DLAP is measured as the mass of analyte actually introduced into the chromatographic column. Ten analytical standards whose concentrations were equally spaced from 0 to 0.45 $\mu\text{g/mL}$ were prepared. The standard containing 0.45 $\mu\text{g/mL}$ represented approximately 10 times the baseline noise for both analytes. These solutions were analyzed with the recommended analytical parameters (10 μL injection). The data obtained were used to determine the required parameters (A and SEE) for the calculation of the DLAP. These parameters and the calculated DLAP's are listed below.

Table 4.2.1
Summary of the Calculated A , SEE , and DLAP

	maneab	zineb
A (ng^{-1})	292800	244900
SEE	82145.4	69878
DLAP (ng)	0.84	0.86

Table 4.2.2
Detection Limit of the Analytical Procedure
for maneb

concentration ($\mu\text{g/mL}$)	mass on column (ng)	peak area
0.000	0.00	0
0.044	0.44	0
0.089	0.89	223710
0.133	1.33	418238
0.177	1.77	487615
0.221	2.21	639389
0.266	2.66	710103
0.310	3.10	871326
0.354	3.54	958243
0.398	3.98	977354
0.443	4.43	1402564

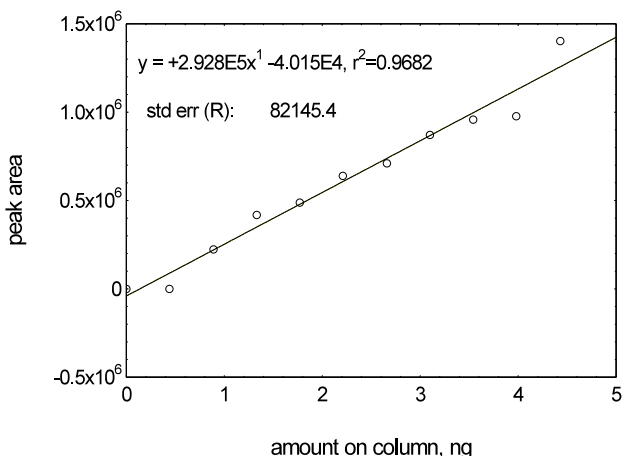


Figure 4.2.1. Plot of the data for determining the DLAP of maneb.

Table 4.2.3
Detection Limit of the Analytical Procedure
for zineb

concentration ($\mu\text{g/mL}$)	mass on column (ng)	peak area
0.000	0.00	0
0.045	0.45	229343
0.090	0.90	320101
0.135	1.35	567710
0.180	1.80	476021
0.225	2.25	560070
0.270	2.70	782160
0.315	3.15	904095
0.360	3.60	956743
0.405	4.05	1081438
0.450	4.53	1165261

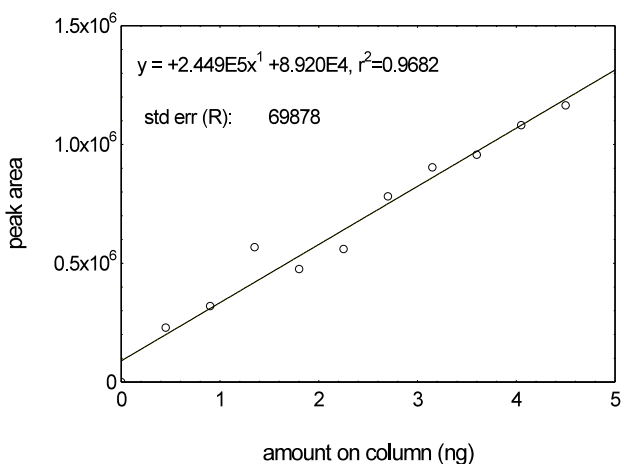


Figure 4.2.2. Plot of the data used for determining the DLAP of zineb.

4.3 Detection limit of the overall procedure (DLOP)

The DLOP is measured as mass per sample and expressed as equivalent air concentration, based on the recommended sampling parameters. Ten samples of maneb/sucrose mixture, ranging in weight from 0.23 to 1.45 mg, were weighed in glass vials containing a mixed-cellulose membrane filter. Ten samples of zineb/sucrose mixture, ranging in weight from 0.17 to 1.22 mg, were similarly prepared. The latter amount, when spiked on a sampler, would produce a peak approximately 10 times the baseline noise for a sample blank. These samples were analyzed with the recommended analytical parameters, and the data obtained used to calculate the required parameters (A and SEE) for the calculation of the DLOP. The parameters obtained and the calculated DLOP's for maneb and zineb are listed below.

Table 4.3.1
Summary of the Calculated A, SEE, and DLOP

	maneb	zineb
A (mg^{-1})	6188000	7635000
SEE	679058	355819
DLOP (mg, sucrose mixture)	0.33	0.14
DLOP (μg)	33	15

Table 4.3.2
Detection Limit of the Overall Procedure
for maneb/sucrose mixture

mass per sample (mg)	peak area
0.23	1909944
0.27	2352685
0.44	3947111
0.47	4418639
0.48	4357205
0.65	4872006
0.94	8026909
0.96	6433869
1.38	8594458
1.45	10464932

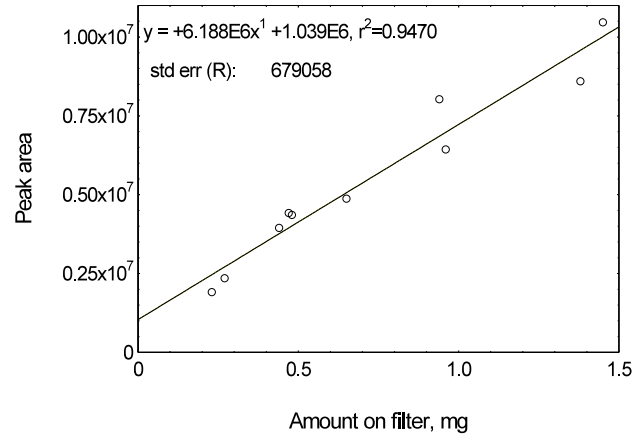


Figure 4.3.1. Plot of data used to determine the DLOP and RQL of maneb/sucrose mixture.

Table 4.3.3
Detection Limit of the Overall Procedure
for zineb/sucrose mixture

mass per sample (mg)	peak area
0.17	2120284
0.18	2057153
0.46	4588090
0.49	3962535
0.50	4285499
0.62	4826301
0.83	6896308
0.94	8300629
1.20	9614353
1.22	10055146

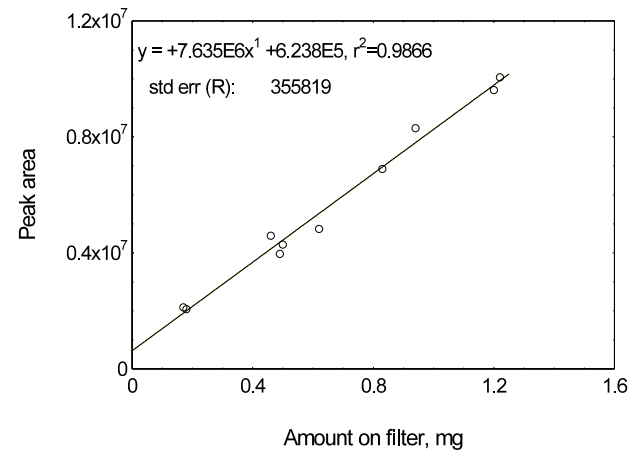


Figure 4.3.2. Plot of data used to determine the DLOP and RQL of zineb/sucrose mixture.

4.4 Reliable quantitation limit

The RQL is considered the lower limit for precise quantitative measurements. It is determined from the regression line data obtained for the calculation of the DLOP (Section 4.3), providing at least 75% of the analyte is recovered. The RQL is defined as the amount of analyte that gives a response (Y_{RQL}) such that

$$Y_{RQL} - Y_{BR} = 10(SD_{BR})$$

therefore

$$RQL = \frac{10(SEE)}{A}$$

The calculated RQL's for maneb and zineb, together with the recoveries at these levels, are listed below. The recoveries are above 75%.

Table 4.4.1
Summary of the RQL's and the Recoveries

	maneb	zineb
RQL (mg, sucrose mixture)	1.10	0.47
RQL (mg/sample)	0.110	0.051
RQL ($\mu\text{g}/\text{m}^3$)	220	103
Recovery (%)	110.1	89.5

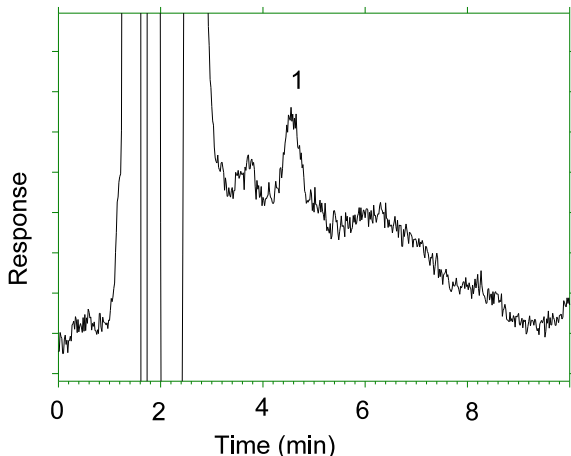


Figure 4.4.1. Chromatogram of the RQL for maneb.
Key: 1 = ethylenebisdithiocarbamate.

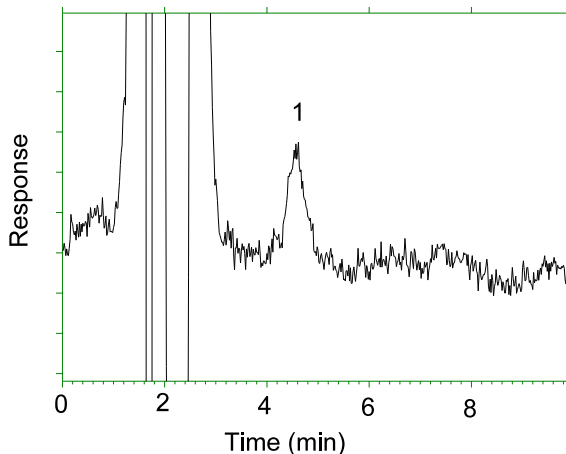


Figure 4.4.2. Chromatogram of the RQL for zineb.
Key: 1 = ethylenebisdithiocarbamate.

4.5 Precision (analytical method)

The precision of the analytical procedure is defined as the pooled relative standard deviation (RSD_p). Relative standard deviations were determined from six replicate injections of analytical standards at 0.5, 0.75, 1, 1.5, and 2 times the target concentration. After assuring that the RSDs satisfy the Cochran test for homogeneity at the 95% confidence level, RSD_p was calculated.

Table 4.5.1
Instrument Response to Maneb

× target concn µg/mL	0.5 × 312	0.75 × 468	1 × 624	1.5 × 936	2 × 1248
peak area	4371777	6643953	8844201	13407652	17642445
	4324183	6578654	8825217	13334192	17818428
	4347347	6572769	8795840	13405380	17738352
	4395817	6676754	8850448	13380629	17884526
	4339860	6694122	8762309	13357394	18026538
	4368155	6622952	8695560	13279650	17753190
avg	4357857	6631534	8795596	13360816	17810580
std dev	25726	49879	58930	48745	133377
RSD %	0.59	0.75	0.67	0.36	0.75

Table 4.5.2
Instrument Response to Zineb

× target concn µg/mL	0.5 × 308	0.75 × 462	1 × 616	1.5 × 924	2 × 1234
peak area	4727680	7073563	9507486	14260604	19078820
	4716993	7090721	9517335	14116518	18976644
	4714717	7133576	9486356	14279337	19006479
	4676611	7114085	9522567	14173887	18964245
	4714064	7094464	9529695	14212995	18900652
	4713370	7114914	9553841	14251726	19009588
avg	4710573	7103554	9519547	14215845	18989405
std dev	17461	21392	22540	61574	58935
RSD %	0.37	0.30	0.24	0.43	0.31

The Cochran test for homogeneity requires the calculation of the g statistics according to the following formula:

$$g = \frac{\text{largest } RSD^2}{RSD_{0.5x}^2 + RSD_{0.75x}^2 + RSD_{1x}^2 + RSD_{1.5x}^2 + RSD_{2x}^2}$$

The *g* statistics obtained were: 0.2741 and 0.3269 for maneb and zineb, respectively. Since these *g* statistics do not exceed the critical value of 0.5065, the RSDs within each level can be considered equal and they can be pooled (RSD_p) to give an estimated RSD for the concentration range studied.

$$RSD_p = \sqrt{\frac{5(RSD_{0.5x}^2 + RSD_{0.75x}^2 + RSD_{1x}^2 + RSD_{1.5x}^2 + RSD_{2x}^2)}{5+5+5+5+5}}$$

The pooled relative standard deviations are 0.64% and 0.34% for maneb and zineb, respectively.

4.6 Precision (overall procedure)

The precision of the overall procedure is determined from the storage data in Section 4.7. The determination of the standard error of estimate (SEE_R) for a regression line plotted through the graphed storage data allows the inclusion of storage time as one of the factors affecting overall precision. The SEE_R is similar to the standard deviation, except it is a measure of dispersion of data about a regression line instead of about a mean. It is determined with the following equation:

$$SEE_R = \sqrt{\frac{\sum (Y_{obs} - Y_{est})^2}{n - k}}$$

n = total no. of data points
k = 2 for linear regression
k = 3 for quadratic regression
 Y_{obs} = observed % recovery at a given time
 Y_{est} = estimated % recovery from the regression line at the same given time

An additional 5% for pump error (SP) is added to the SEE_R by the addition of variances to obtain the total standard error of estimate.

$$SEE = \sqrt{(SEE_R)^2 + (SP)^2}$$

The precision at the 95% confidence level is obtained by multiplying the standard error of estimate (with pump error included) by 1.96 (the *z*-statistic from the standard normal distribution at the 95% confidence level). The 95% confidence intervals are drawn about their respective regression lines in the storage graphs, as shown in Figures 4.7.1.1 to 4.7.2.2. The precisions of the overall procedure are ±19.8% and ±19.4% for maneb and zineb, respectively.

4.7 Storage test

Storage samples were prepared by weighing the maneb/sucrose (at a level of around 27 mg, or maneb content of 2.5 mg) or zineb/sucrose mixture (at a level of around 24 mg or zineb content of 2.5 mg) in a vial. Thirty-six samples were prepared for each analyte. Six samples were analyzed on the day of preparation. The rest of the samples were divided into two groups: 15 were stored at 5°C, and the other 15 were stored at ambient temperature (about 22°C) in a closed drawer. At 2-4 day intervals, three samples were selected from each of the two storage sets and analyzed.

Table 4.7.1
Storage Test for Maneb

time (days)	percent recovery (ambient)			percent recovery (refrigerated)		
	0	103.0	102.3	96.2	103.0	102.3
0	98.8	99.6	100.1	98.8	99.6	100.1
4	118.3	103.0	114.1	100.6	99.3	103.2
6	102.6	102.7	101.5	91.3	89.8	93.0
8	87.2	99.6	70.8*	81.2	94.1	106.6
12	82.3	77.9	95.5	70.8	96.9	57.7*
15	91.1	103.1	101.4	89.9	91.6	85.9

* Outliers, not used

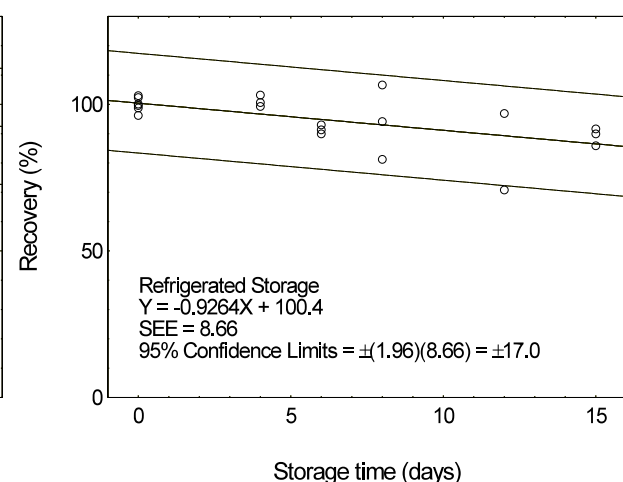
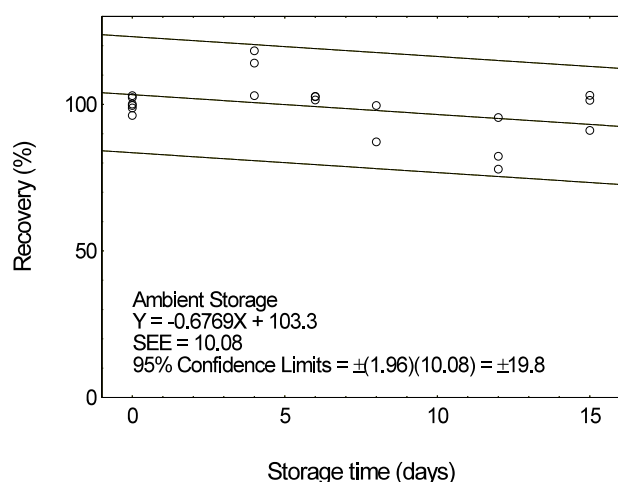


Figure 4.7.1.1. Ambient storage test for maneb.

Figure 4.7.1.2. Refrigerated storage test for maneb.

Table 4.7.2
Storage Test for Zineb

time (days)	percent recovery (ambient)			percent recovery (refrigerated)		
	0	102.1	102.2	103.5	102.1	102.2
0	102.7	103.7	85.9	102.7	103.7	85.9
4	102.1	100.8	98.0	103.2	100.7	89.4
6	90.0	92.5	83.1	90.7	91.5	87.3
8	129.9*	117.4	92.9	77.4	104.4	79.9
12	90.5	104.4	109.4	76.0	101.6	98.4
15	88.9	99.1	93.2	94.6	96.8	94.6

*outlier, not used

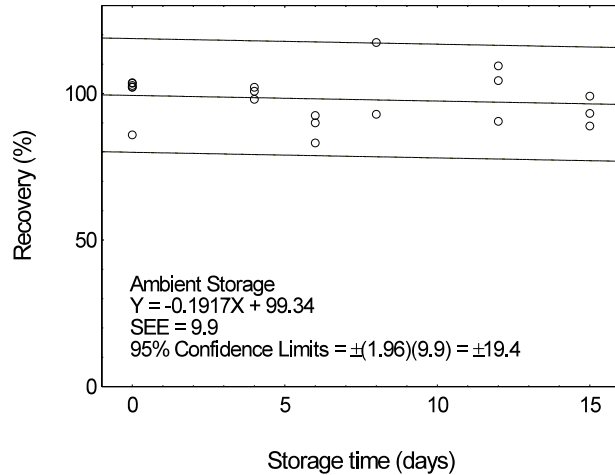


Figure 4.7.2.1. Ambient storage test for zineb.

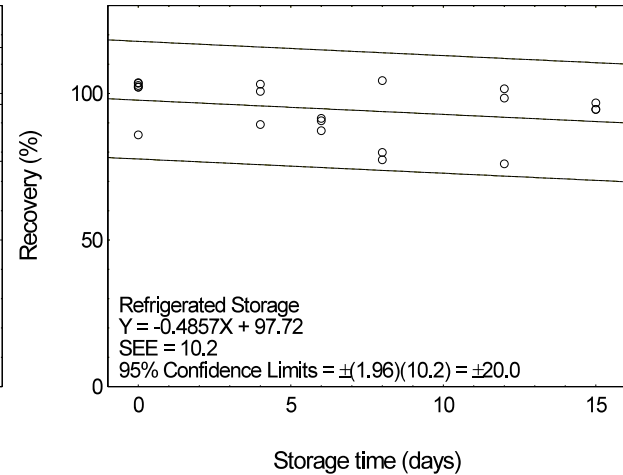


Figure 4.7.2.2. Refrigerated storage test for zineb.

4.8 Reproducibility

Reproducibility samples were prepared by weighing maneb or zineb/sucrose mixtures in a scintillation vials containing a GN-4 filter (membrane filter made from mixed esters of cellulose). The samples were submitted to an SLTC service branch for analysis. The samples were analyzed after being stored for 2 days at ambient temperature. No sample result had a deviation greater than the precisions of the overall procedure determined in Section 4.7, which are $\pm 19.8\%$ and $\pm 19.4\%$ for maneb and zineb, respectively.

Table 4.8.1
 Reproducibility Data for Maneb

mg expected	mg found	percent found	percent deviation
1.355	1.34	98.9	-1.1
1.150	1.05	91.3	-8.7
1.755	1.74	99.1	-0.9
2.370	1.97	83.1	-16.9
1.564	1.47	94.0	-6.0
2.494	2.28	91.4	-8.6

Table 4.8.2
 Reproducibility Data for Zineb

mg expected	mg found	percent found	percent deviation
2.627	2.70	102.8	+2.8
1.310	1.25	95.4	-4.6
1.935	1.66	85.8	-14.2
2.193	2.21	100.8	+0.8
1.827	1.48	81.0	-19.0
1.535	1.41	91.8	-8.2

4.9 Sampler capacity

Sampling capacity was not tested. Maneb and zineb have very low vapor pressure and are not expected to vaporize or sublime significantly at ambient temperature. Generally one would not expect dust particles to break through a membrane filter. Retention efficiencies were tested by pulling 500 L of 80%-RH air through cassettes containing about 12 to 18 mg of maneb/sucrose or zineb/sucrose mixture (containing about 1.2 to 1.8 mg of maneb or zineb). The recoveries were 97.0% and 102.8% for maneb and zineb, respectively.

Table 4.9.1
Retention Efficiency of Maneb

Treated with 500 L of humid air			Control group		
mg spiked*	peak area	area/mg	mg spiked*	peak area	area/mg
13.27	4714261	355257	13.43	4865505	362286
17.87	6473108	362233	16.26	5520337	339504
12.82	4077903	318089	13.94	5024973	360472
15.68	5545443	353663	13.45	4654504	346060
13.56	4833093	356423	17.52	Lost	---
16.37	5162170	315343	14.90	5409830	363076
	average =	343501		average =	354280
Retention efficiency = 343501/354280 = 97.0%					

* Amount of mane/sucrose mixture. Maneb content = 10.1%, by weight.

Table 4.9.2
Retention Efficiency of Zineb

Treated with 500 L of humid air			Control group		
mg spiked*	peak area	area/mg	mg spiked*	peak area	area/mg
16.32	5137892	314822	17.04	5905005	346538
17.08	6340746	371238	17.96	6559759	365243
17.83	6824216	382738	16.76	6171506	368228
17.23	6558921	380669	16.73	6021564	359926
17.67	7121568	403032	16.69	5983878	358531
18.01	7100030	394782	16.24	6304134	388186
	average =	374547		average =	364442
Retention efficiency = 374547/364442 = 102.8%					

* Amount of zineb/sucrose mixture. Zineb content = 11.0%.

4.10 Extraction efficiency and stability of extracted samples

4.10.1 Extraction efficiency

Samples for the extraction efficiencies (EE) of mane and zineb were prepared by weighing, at 0.05 to 2 times the target concentrations, the mane or zineb-sucrose mixture in a scintillation vial containing a GN-4 filter. These samples were stored overnight at ambient temperature and then extracted and analyzed. The average extraction efficiencies over the working range of 0.5 to 2 times the target concentration were 100.4% and 98.5% for mane and zineb, respectively.

Table 4.10.1.1
Extraction Efficiency (%) for Maneb

0.05×TC		0.1×TC		0.2×TC		0.5×TC		1×TC		2×TC	
mg*	EE	mg*	EE	mg*	EE	mg*	EE	mg*	EE	mg*	EE
1.45	89.7	2.56	90.2	5.26	112.0	13.11	109.2	25.61	104.1	51.67	100.4
1.30	122.0	2.29	100.9	5.23	79.9	12.77	102.0	24.48	105.5	50.43	102.8
1.29	83.6	2.27	104.1	5.13	65.2	12.71	96.3	25.46	110.1	50.21	97.6
1.21	92.4	2.19	87.3	4.72	95.8	12.13	95.4	24.96	98.4	49.99	97.9
1.13	96.7	2.14	115.1	4.55	94.4	11.75	97.2	23.73	104.9	48.26	95.5
1.07	97.2	2.05	106.0	4.29	106.8	11.62	97.7	23.61	94.4	47.69	99.0
\bar{X}	96.9	\bar{X}	100.6	\bar{X}	92.3	\bar{X}	99.6	\bar{X}	102.9	\bar{X}	98.9

* Amount of mane/sucrose mixture. Maneb content = 10.1%.

Table 4.10.1.2
Extraction Efficiency (%) for Zineb

0.05×TC		0.1×TC		0.2×TC		0.5×TC		1×TC		2×TC	
mg*	EE	mg*	EE	mg*	EE	mg*	EE	mg*	EE	mg*	EE
1.30	88.2	2.73	82.3	5.47	97.1	13.03	100.3	25.60	98.9	51.16	100.3
1.29	87.7	2.69	104.3	5.11	99.0	12.79	101.8	25.35	96.0	50.35	94.9
1.29	106.8	2.64	108.2	5.01	110.6	12.67	101.6	25.17	99.2	50.12	99.3
1.21	98.1	2.33	95.8	4.82	102.4	12.47	96.1	24.43	92.9	49.79	99.8
1.18	101.0	2.25	117.4	4.57	88.6	12.46	99.5	24.31	93.8	49.51	101.9
1.05	97.0	2.22	108.5	4.52	106.0	12.43	98.3	24.24	97.9	49.08	100.3
\bar{X}	96.5	\bar{X}	102.7	\bar{X}	100.6	\bar{X}	99.6	\bar{X}	96.4	\bar{X}	99.4

* Amount of zineb/sucrose mixture. Zineb content = 11.0%.

4.10.2 Stability of extracted samples

The stability of the extracted samples was investigated by reanalyzing the target concentration samples 24 h after initial analysis. After the original analysis was performed, three vials were recapped with new septa while the remaining three retained their punctured septa. The samples were reanalyzed with fresh standards.

Table 4.10.2.1
Stability of extracted samples for maneb

punctured septa replaced			punctured septa retained		
initial EE (%)	EE after one day (%)	difference	initial EE (%)	EE after one day (%)	difference
104.1	103.0	-1.1	98.4	98.5	+0.1
105.5	104.0	-1.5	104.9	103.8	-1.1
110.1	111.6	+1.5	94.4	97.8	+3.4
averages			averages		
106.6	106.2	-0.4	99.2	100.0	+0.8

Table 4.10.2.2
Stability of extracted samples for zineb

punctured septa replaced			punctured septa retained		
initial EE (%)	EE after one day (%)	difference	initial EE (%)	EE after one day (%)	difference
98.9	97.7	-1.2	92.9	92.5	-0.4
96.0	94.8	-1.2	93.8	94.2	+0.4
99.2	99.4	+0.2	97.9	96.5	-1.4
averages			averages		
98.0	97.3	-0.7	94.9	94.4	-0.5

4.11 Qualitative analysis

As an alternative analytical procedure, samples of maneb or zineb can be analyzed by atomic absorption spectrometry (AA), if the samples are well digested in acid. The chelated metals do not give full-strength response on AA. The Inductively Coupled Plasma (ICP) emission spectrometry may also be used.

Afsar and Demirata (Ref. 5.12) reported a method of differentiating maneb, zineb, and mancozeb on the basis of colors produced after treatment of saturated solutions of the fungicides in propanol-acetone (1:1, v/v), first with dithizone and then with monobasic sodium phosphate solution in the same solvent.

One can also analyze by gas chromatography the released carbon disulfide after the acid-decomposition of maneb or zineb. However, there is a danger of interference from thiram, a related fungicide that is sometimes used together with maneb or zineb.

4.12 Stability of the EBDTC in the extraction solvent

Three sets of stock solutions were prepared in the extraction solvents of various freshness: 1-, 27-, and 49-day old. Their instrument responses were followed periodically for 5 days after the preparation. The results are plotted in the following graphs. The instrument responses remained essentially constant during this period.

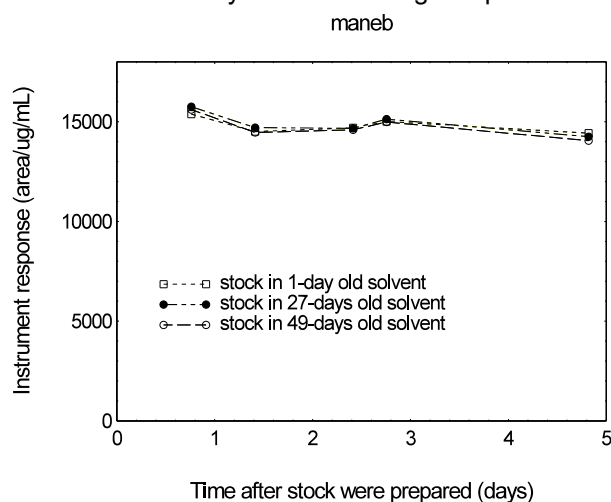


Figure 4.12.1. Stability of maneb stock solutions prepared in the extraction solvents of various age.

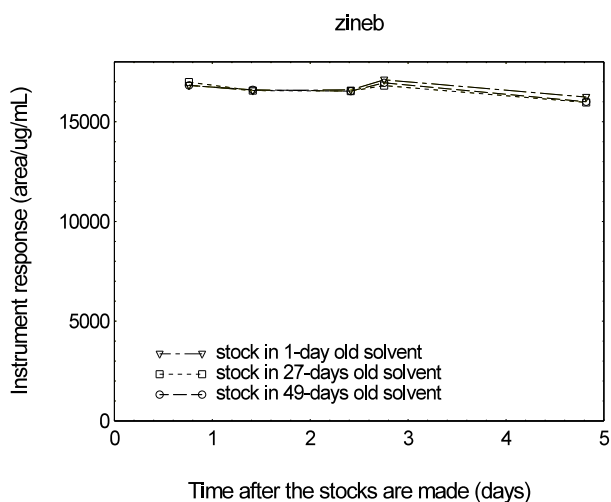


Figure 4.12.2. Stability of zineb stock solutions prepared in the extraction solvents of various age.

5. References

- 5.1. *Thiocarbamate Pesticides: A General Introduction, Environmental Health Criteria 76*, World Health Organization, Geneva, 1988.
- 5.2. *IARC Monographs on the Evaluation of Carcinogenic Risk of Chemicals to Man*. Geneva: World Health Organization, International Agency for Research on Cancer, 1972-1985.
- 5.3. Newsome, W.H., "Ethylenebisdithiocarbamates and Their Degradation Products", in *Analytical Methods for Pesticides and Plant Growth Regulators*, Vol. XI, Academic Press, NY, 1980.
- 5.4. *Merck Index*, Windholz, M., Ed., 10th ed., Merck & Co., Rahway, NJ, 1983.
- 5.5. *Merck Index*, Budavari, S., Ed., 11th ed., Merck & Co., Rahway, NJ, 1989.
- 5.6. Bardarov, V. and C. Zaikov, *J. Chromatogr.*, **1989**, 479, 97-105.
- 5.7. Miles, C.J. and M. Zhou, *J. Assoc. Off. Anal. Chem.*, **1991**, 74(2), 384-8.
- 5.8. Kunugi, Y., *Shiryō Kenkyū Hokoku (Tokyo Hishiryō Kensasho)*, **1994**, 19, 56-57.
- 5.9. Jongen, M.J.M., J.S. Ravensberg, R. Engel, and L.H. Leenheers, *J. Chromatogr. Sci.*, **1991**, 29(7), 292-7.
- 5.10. Marple, V.A., B.Y.H. Liu, and K.L. Rubow, *Am. Ind. Hyg. Assoc. J.*, **1978**, 39, 26-32.
- 5.11. Gosselin, R.E., R.P. Smith, H.C. Hodge. *Clinical Toxicology of Commercial Products.*, 5th ed. Baltimore: Williams and Wilkins, 1984.

5.12. Afsar, H. And B Demirata, *J. Assoc. Off. Anal. Chem.*, **1987**, 70(5), 923-4.